

QUICKIE 2

ONE SIZE DOES NOT FIT ALL.
BUT ONE CHAIR DOES.



THE FOLDING CHAIR THAT FITS JUST RIGHT

The most popular ultra lightweight folding wheelchair—engineered to adapt over time. With over 100,000 configuration options, it delivers personalized performance for a wide range of riders and needs.

TRUE FIT

Includes one free growth kit within five years, making it ideal for pediatrics or anyone needing flexibility. Frame sizes start at 12"x12", with a depth-adjustable rear frame for added adaptability.

LIGHT, STRONG, RELIABLE

A diamond-shaped crossbrace delivers added rigidity and long-term durability. Crafted from 7000 series aerospace aluminum—the same material used in rockets—this unique design offers exceptional strength without added weight.

FITS TODAY AND TOMORROW


Offered in 11 seat depths, 13 frame widths, and 8 frame depths. Caster housings, forks, and wheel sizes can be tailored to optimize performance across various environments.

MEETS SIMPLE TO COMPLEX SEATING AND MOBILITY NEEDS

Supports a wide range of seating, from basic comfort to advanced clinical positioning. Pair with JAY cushions and backs, including air, fluid, and cooling options, for personalized pressure relief and support.

POWER WHEN NEEDED

The strong, durable frame supports a stable ride and is compatible with all Empulse power assist devices. R90 aids independent propulsion, R20 supports caregivers over distance or hills, and F35 adds powered pull in front for varied terrain.

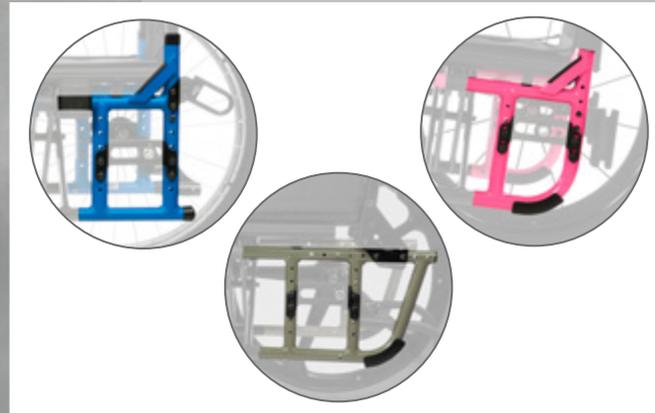

TECHNICAL DATA

HCPCS CODE	K0005
FRAME MATERIAL	7000 Aluminum
SEAT WIDTH	12" to 24"
SEAT DEPTH	12" to 22"
FRAME ANGLES	60°, 70°, 80°, 90°
SEAT HEIGHT FRONT	13.5" to 22"
SEAT HEIGHT REAR	13" to 22"
BACK HEIGHT ADJUSTMENT	9" to 21"
BACK ANGLE ADJUSTMENT	From 87° to 105°
CENTER OF GRAVITY (COG)	1" to 3.5"
CHAIR WEIGHT	25.4 lbs.
USABLE TRANSPORT WEIGHT	17.8 lbs.
WEIGHT CAPACITY	Standard: 300 lbs., Heavy Duty: 350 lbs.
TRANSIT	Yes WC19 Transit approved.

Refer to order form as specifications may vary based on configuration.

Learn more



This document can be viewed in PDF format at www.SunriseMedical.com.



Sunrise Medical (US)
MK-100079 Rev.F ©05.2025



800.333.4000

www.SunriseMedical.com

