


zippie<sup>®</sup>  
iris

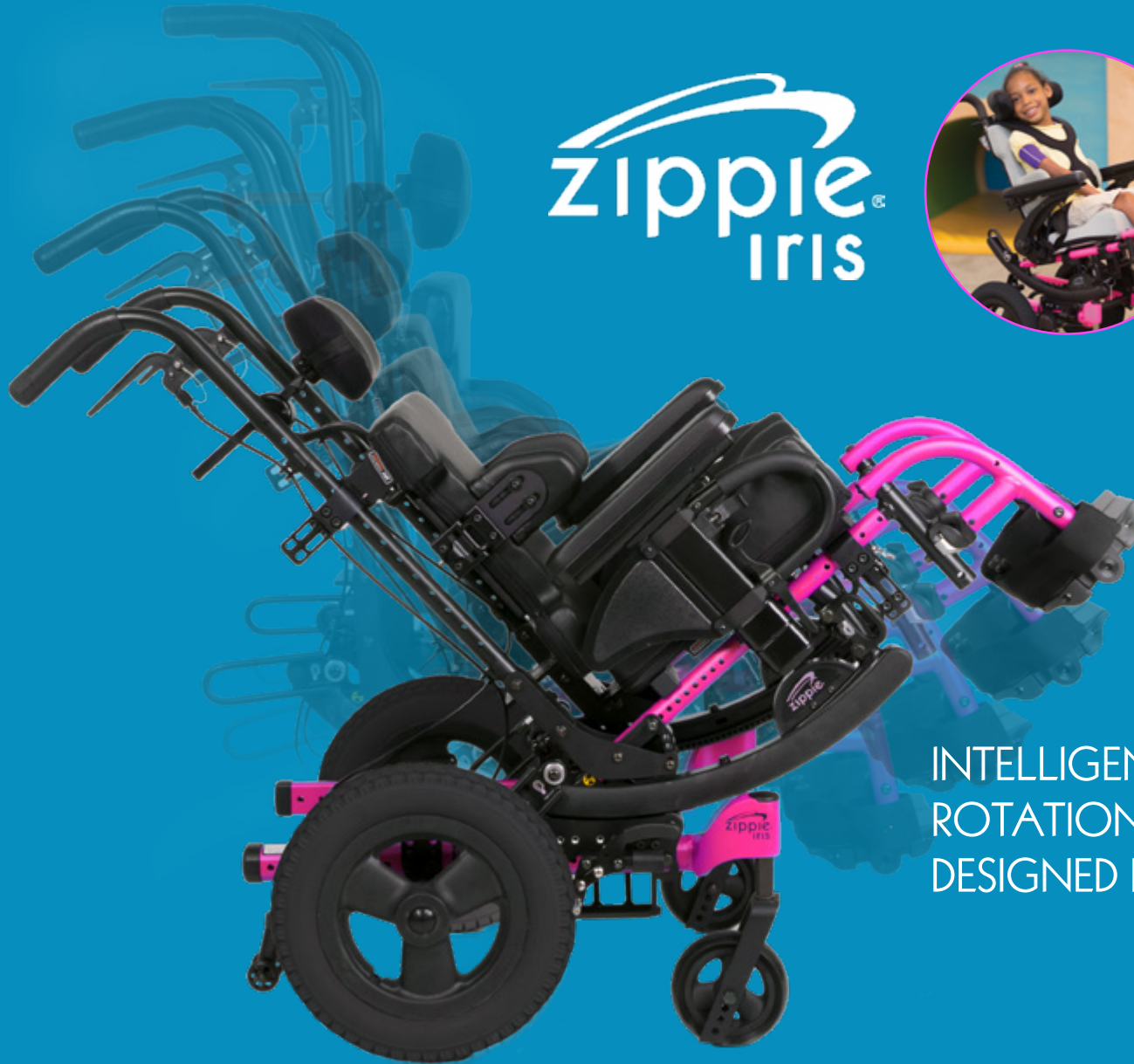

INTELLIGENT  
ROTATION-IN-SPACE  
DESIGNED FOR KIDS

## 40° OR 55° OF INTELLIGENT ROTATION IN SPACE TECHNOLOGY

The ZIPPIE IRIS® features patented Intelligent Rotation in Space technology, allowing the seat frame to rotate around the child's center of gravity. Designed specifically for pediatric riders, this ensures smooth, secure tilting and easy operation for caregivers. Choose from 40° or 55° rotation ranges, adjustable to provide up to 45° or 60° posterior tilt for superior positioning.

## GROWTH CAPABILITIES

Designed to grow with your child, the IRIS folding frame includes growth up to 4" of depth and 2" of width. With the True Fit Growth program, the original owner receives one free growth kit within five years (not available with SE models).

## Z-FINITY FOOTREST SYSTEM

Supports evolving lower extremity needs, adjusting from 0.5" above to 21" below the seat pan. Its rotating hardware allows infinite footrest angle and placement changes for optimal foot, ankle, and knee positioning.

## OPTIONAL XLOCK® FOLDING CROSSBRACE

Enjoy the benefits of both rotation-in-space and folding convenience. The optional XLOCK crossbrace locks when open for rigidity and folds compactly for transport.

## ACCESSORIES FOR CAREGIVER CONVENIENCE

Caregivers benefit from features like the Foot Release Tilt Actuator (no exposed cables), optional easy-to-use push button hub locks, and an Angle Adjustable Stroller Handle with 260° of adjustability. A height-adjustable backrest comes standard.

## VARIOUS SEATING OPTIONS

Customize the IRIS with dynamic or rigid seating. The Dynamic Back accommodates movement to reduce damage and maintain positioning. The MONO Backrest System offers lateral support, flexibility, and an optional dynamic feature. Compatible LECKEY seating systems, such as the Squiggles or Mygo, provide targeted postural support.



## SPECIFICATIONS

HCPCS CODE	RIGID: E1233 FOLDING: E1234
PRODUCT WIDTH	17.5" to 25.5"
SEAT WIDTH	RIGID 10" to 18", Folding: 10" to 16" (grows to 18")
SEAT DEPTH	10" to 18"
TILT RANGE	40° / 55°
FRONT SEAT-TO-FLOOR HEIGHT	13" to 19"
REAR SEAT-TO-FLOOR HEIGHT	13" to 19"
BACK ANGLE ADJUSTMENT	85° to 120° (5° fore to 30° aft)
TRANSIT	Yes*
RIGID STANDARD WEIGHT	32 lbs., Transport Weight**: 22.8 lbs
FOLDING STANDARD WEIGHT	35.8 lbs., Transport Weight**: 26.6 lbs.
WEIGHT CAPACITY	225 lbs. (165 lbs. with Z-Finity)

\* Conforms with ANSI/RESNA WC-19 standards. Please refer to the QUICKIE & ZIPPIE WC-19 Transit Securement Supplement (MK-100072) for additional information.

A supplement is available

\*\* Excludes footrests, rear wheels, seat pan, armrests, and transit.



Conforms to ANSI/RESNA  
WC-19 transit standards

Learn more



This document can be viewed in PDF  
format at [www.SunriseMedical.com](http://www.SunriseMedical.com).



MK-100293 Rev C ©09.2025



800.333.4000

[www.SunriseMedical.com](http://www.SunriseMedical.com)

