

Manta Classic Back Support

The Spex Manta Classic back support combines the best of two Spex Seating favourites, merging the familiar, wraparound shape of the Manta with the proven cube contouring technology of the Spex Classic. This has opened up the Manta to a wider range of users who may benefit from more customisable and targeted support.

On the inside are built-in lateral trunk supports, a pressure-relieving foam inner and the ability to contour on the spot. On the outside, a soft, breathable cover hugs the whole back support, offering comfort to the user.

CLINICAL

The Manta Classic is designed for individuals who require:


- Contouring to provide trunk alignment and accommodate mild to moderate trunk asymmetries
- Comfortable shoulder-height back support to relieve pressure and reduce pain
- Adjustable lateral trunk supports to enhance trunk stability and spinal alignment
- Flexibility to address evolving postural needs
- Improved upper body function and independence due to minimal hardware interference around lateral supports
- Increased comfort around bony prominences through combined pressure-relieving foam and adjustable cube immersion

FUNCTIONAL

- Adjustable Cube Contouring: Add or remove cubes while the user remains seated. Simply unzip the cover from the top to access the contouring system
- Cube Anchor System: Prevents shifting or slipping of cubes during transfers
- Wing Cube Pockets: Contouring cubes for the lateral wings are included in the box and can be added as needed
- Extended Lateral Hardware: Maintains the Manta's signature shape, even with built-up contouring
- Pressure-Relieving Inner Overlay for maximum user comfort
- Neoprene Lateral Covers: Lateral trunk support covers are made from soft neoprene, ensuring comfort during self-propulsion
- Integrated Harness Retainer (Please note this is not intended for clinical positioning – for this purpose, please add Spex Harness Guides)



SIZE CHART

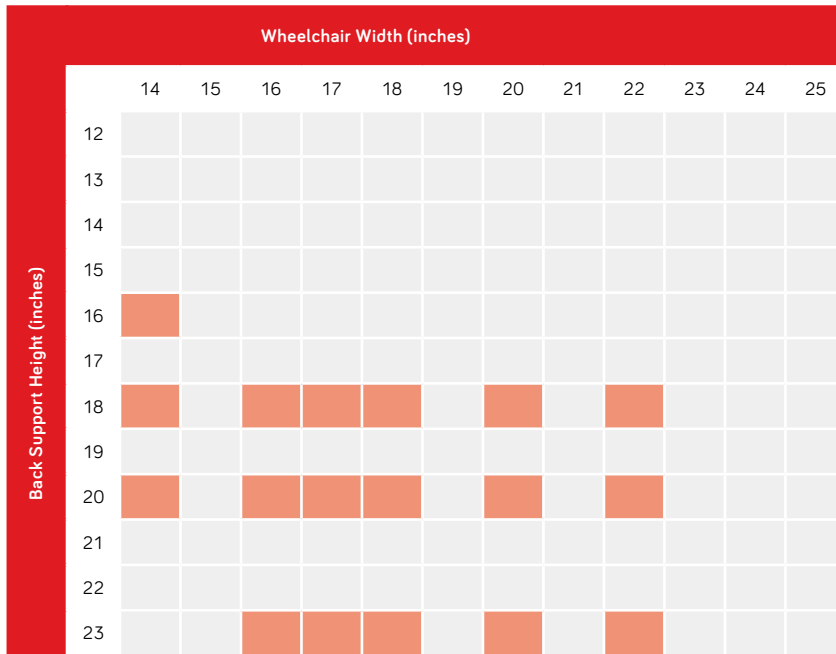

- Red squares are standard sizes.
- Non-standard or custom sizes are available on request.

PRODUCT PERFORMANCE

Postural Control	5
Lateral Stability	5
Postural Asymmetries	5
Changing Needs	4
Pressure Care	4

• Note: This is a guide only and ratings are dependent on correct use of product.

PRODUCT DETAILS

Width	Capacity	Product Weight
10"-15"	200lbs (90kg)	(16x16)
16"-20"	300lbs (136kg)	2.2kg

• Note: This is the weight of the back support only and does not include the hardware. For hardware weight, refer to p28.

MANTA FITTING INFORMATION

Back Support Width	Min. Internal Width (A)	Max. Internal Width (B)	Lateral Wing Size	Lateral Height Adjustment	Lateral Depth Adjustment
10"	5.5"	6.5"	0	3cm	5cm
12"	7.5"	8.5"			
13"	8"	11"			
14"	9"	12"	1	3cm	5cm
15"	10"	13"			
16"	11"	14.5"			
17"	12"	15.5"	2	4cm	7.8cm
18"	13"	16.5"			
19"	14"	17.5"			
20"	15"	18.5"			
22"	17"	20.5"			

