


GATEWAY[®] 3G



Solid Surface Portable Ramp

Designed for safety and engineered for reliability, the GATEWAY 3G, a solid surface portable wheelchair ramp, is a smart alternative to costly wooden ramps. When mobility needs extend beyond the scope of a conventional portable ramp, this durable ramp is your answer. Available with or without handrails.

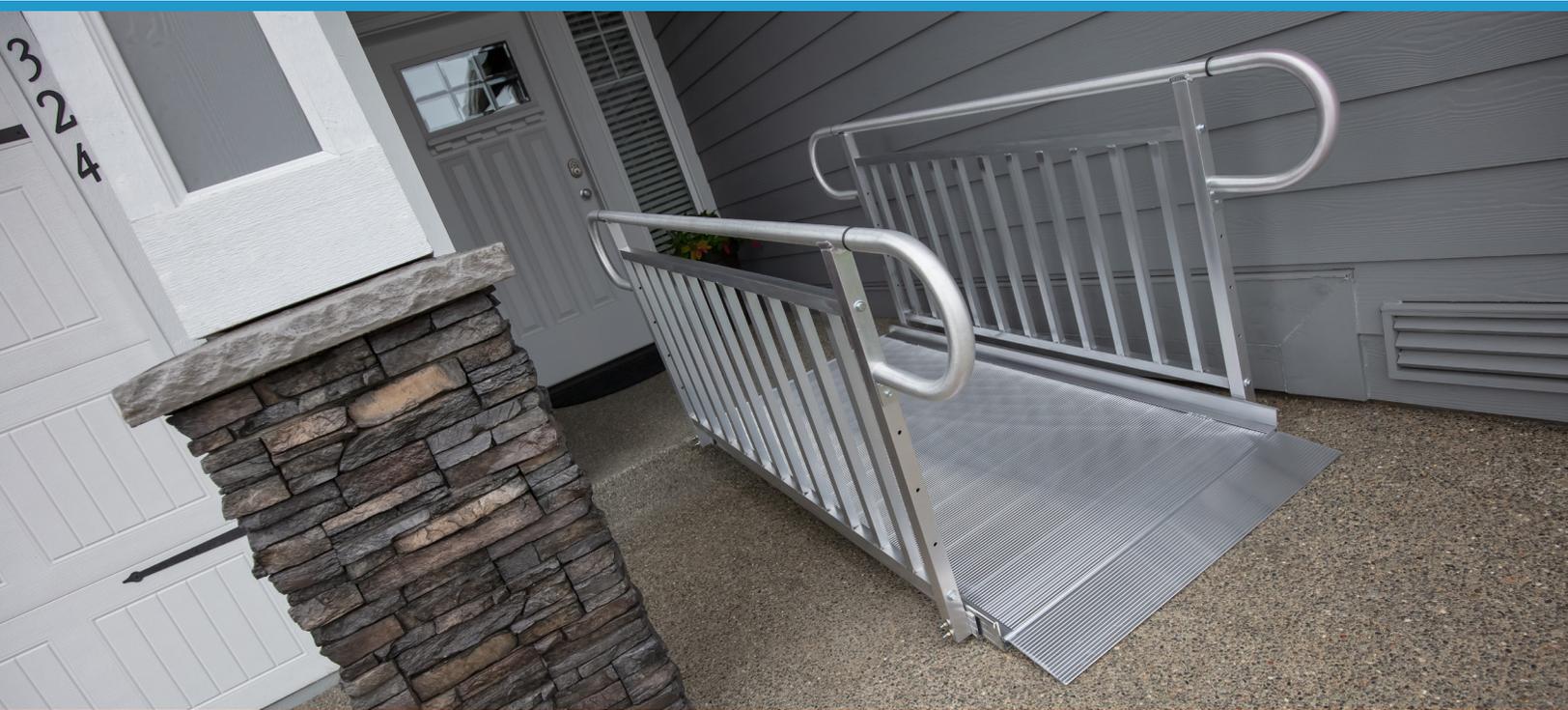

The GATEWAY 3G Advantage

- A safe ramp solution designed to provide home access. Available with or without handrails.
- Can be used for semi-permanent applications and can easily be moved or relocated as needs change.
- With lengths up to 10', the GATEWAY 3G can accommodate a series of entryway steps.
- Features a formed top transition for secure placement, and an adjustable bottom transition that allows smooth ramp-to-ground conversion.
- Designed for the rigors of repeated, everyday use in all outdoor weather conditions but can also be used indoors.
- Durable construction for years of maintenance-free use. Will not rust, rot, or deteriorate.
- The solid, seam-free surface features permanent (extruded) tread for slip-resistance and superb traction.
- Raised side rails, also known as curbs, help to prevent mobility devices from driving off the edge of the ramp.
- Made of lightweight, heavy-duty aluminum and features a 1,000 lb. weight capacity.
- Available in 3-, 4-, 5-, 6-, 7-, 8-, 9-, & 10-foot lengths with a standard width of 36 inches.

Visit ezaccess.com for product specifications

Safety is our top priority, and we are committed to making the process of choosing the best ramp solution as worry-free as possible. With EZ-ACCESS, you can be confident that your accessibility needs will be resolved safely and quickly.



Ramp Only

- A great option for users who don't require handrail assistance when making their way up or down the ramp.
- The 1-piece design means no assembly or movable parts; just secure the ramp in the desired position, and it's ready to be used!



Ramp with Two-line Handrails and Loops

- Two-line handrails with loops offer increased safety and user stability.
- Features smooth, continuous top rails - which are welded for extra strength and durability - and a lower second rail offers an additional gripping surface.
- Assembly is quick and easy, requiring only minimal tools to attach the handrail posts and loops.



Ramp with Vertical Pickets and Loops

- Smooth, continuous top rails with vertical pickets and loops offer an added level of safety, especially for children and pets. Pickets meet the 4-inch sphere rule for spacing.
- Picketed insert sections are pre-assembled for quick and secure installation.
- Assembly is quick and easy, requiring only minimal tools to attach the handrail posts and loops.

All GATEWAY® 3G Solid Surface Portable Ramps are covered under the EZ-ACCESS lifetime warranty and are proudly made in the USA.



customerservice@ezaccess.com | 1.800.451.1903 | www.ezaccess.com